


Reinforced PVC Endotracheal Tube With Suction Lumen

Reinforced Endotracheal Tube with Suction Lumen is a single-use airway management solution made of medical-grade PVC¹ with excellent chemical stability.

- The airway tube is smooth inside and outside, with no burrs and no impurities. Latex and DEHP free.
- The Bevel tip and Murphy hole are smoothed to minimize damage to the airway during insertion.
- X5CrNi18-10 austenitic stainless steel of spring coil, kink-free tube eliminates the risk of airway tube occlusion.
- X-ray impermeable blue line design for easy identification of tube tip position.
- The suction lumen effectively aspirates subglottic secretions and reduces the incidence of iatrogenic infections and ventilator-associated pneumonia.
- Curvature design of compliant tracheal tube² (Magill Curve).
- Thin-walled, soft and frosted cuff, with a smooth transition of connecting surfaces and a high-capacity, low-pressure design that reduces damage to the airway.

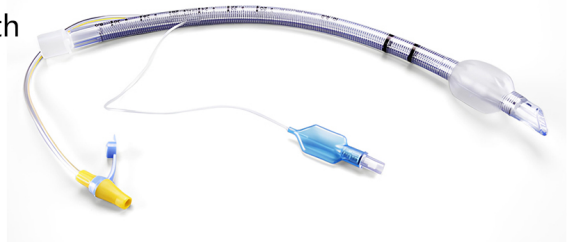

Pro manufacturer of high-end medical consumables

Medical grade UV adhesive, excellent tensile strength.³

Highly elastic airway tube with excellent flexibility and compression resistance.⁴
(Shore hardness from 20A to 65A can be customized).

The left-facing bevel tip improves the visualization of the vocal cords during laryngoscopy. It also increases the diameter of the primary opening, making it less likely to be blocked by secretions.


1 : Silicone-containing PVC or TPU material available. PU cuff is optional. Cool white

2 : Precise design of the curvature of the airway tube using 3D scanning technology.

3 : Failure testing of individual connection components using tensile and compression testers ensures that the tensile strength meets design criteria.

4 : Bite force simulation of the main tube using a compression tester ensures that the compression resistance meets the design criteria.

Reinforced PVC Endotracheal Tube With Suction Lumen

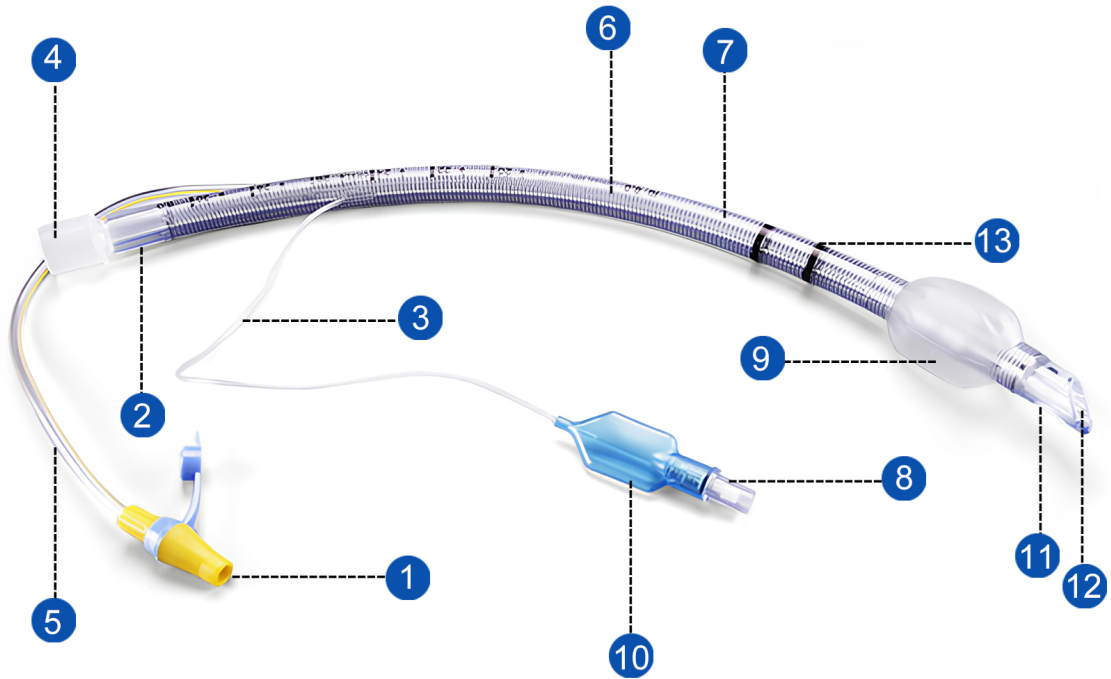

- | | | | |
|----------------------------------|------------------------|-------------------|-------------------|
| 1) Suction Connector | 2) Radio-opaque Line | 3) Inflation Line | 4) 15mm Connector |
| 5) Suction Line | 6) Airway Tube | 7) Spring Coil | 8) Valve |
| 9) High Volume Low Pressure Cuff | 10) Pilot Balloon | 11) Murphy Eye | |
| 12) Bevel & Tip | 13) Depth Marker Lines | | |

Size	5.0#	5.5#	6.0#	6.5#	7.0#	7.5#	8.0#	8.5#	9.0#
Packing:	20pcs/box				10pcs/box				

Widely used in the following applications:

- Intensive Care Unit
- Oral & Maxillofacial Surgery
- Ambulance Equipment
- Emergency Dept.
- Respiratory Dept.
- Anesthesiology Dept.
- Pre-hospital Emergency

Color: Cool white
 Size: Cool white