


## Single Lumen Classic Laryngeal Mask Airway

This classic Laryngeal Mask Airway is a versatile, single-use<sup>1</sup> airway made of biocompatible and non-toxic medical silicone.

- With non-adhesive debridement process, smooth shaft and more comfortable to use. Less irritation and cardiovascular response than tracheal tube.
- The cuff is soft and has a matte finish that reduces the incidence of seal leaks and shifting. It is molded in one piece with a non-adhesive back and can be customized with aperture bars or without.
- Multiple sizes available, suitable for neonate, infant, child and adult.
- Designed to fit the curvature of the human airway.<sup>2</sup>
- Flesh, blue or clear color of cuff is optional

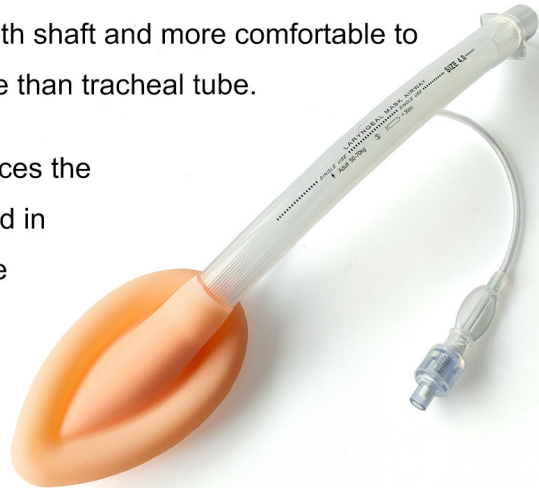

Pro manufacturer of high-end medical consumables

Fit with human oral physiology structure with little tissue damage, easy to insert.

Great compression resistance of the airway tube.<sup>3</sup>

Excellent flexibility of the cuff with superior sealing performance.

Medical grade UV adhesive, excellent tensile strength.<sup>4</sup>

1, Reusable or single-use can be customized.


2, Precise design of the curvature of the airway tube using 3D scanning technology.

3, Bite force simulation of the main tube using a compression tester ensures that the compression resistance meets the design criteria.

4, Failure testing of individual connection components using tensile and compression testers ensures that the tensile strength meets design criteria.



## Classic Laryngeal Mask Airway Specification

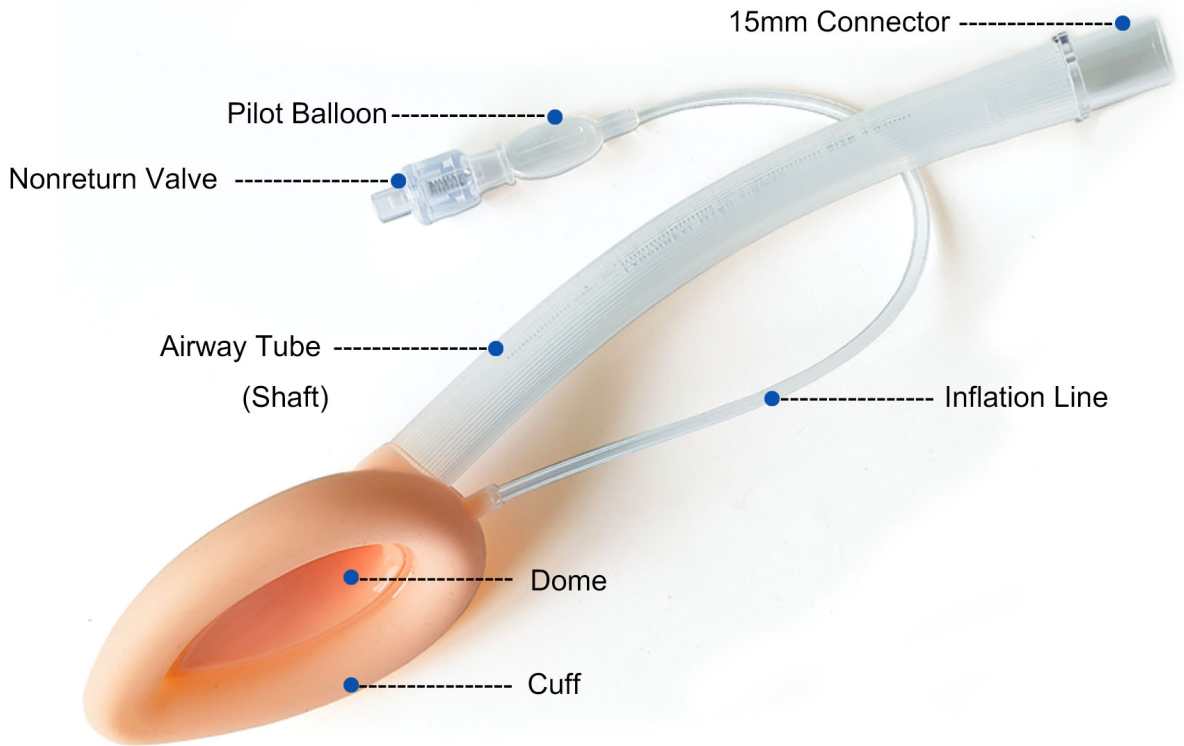

Size	1.0#	1.5#	2.0#	2.5#	3.0#	4.0#	5.0#
Patient Weight	<5kg	5-10kg	10-20kg	20-30kg	30-50kg	50-70kg	70-100kg
Cuff volume	<4ml	<7ml	<10ml	<14ml	<20ml	<30ml	<40ml
Package	100pcs/ctn						

\*It is recommended that the cuff be inflated to a maximum of 60 cm H<sub>2</sub>O intracuff pressure.

Widely used in the following applications:

- Ambulance Equipment
- Emergency Dept.
- Respiratory Dept.
- Anesthesiology Dept.
- Pre-hospital Emergency
- Ambulance Vehicles